



Launch

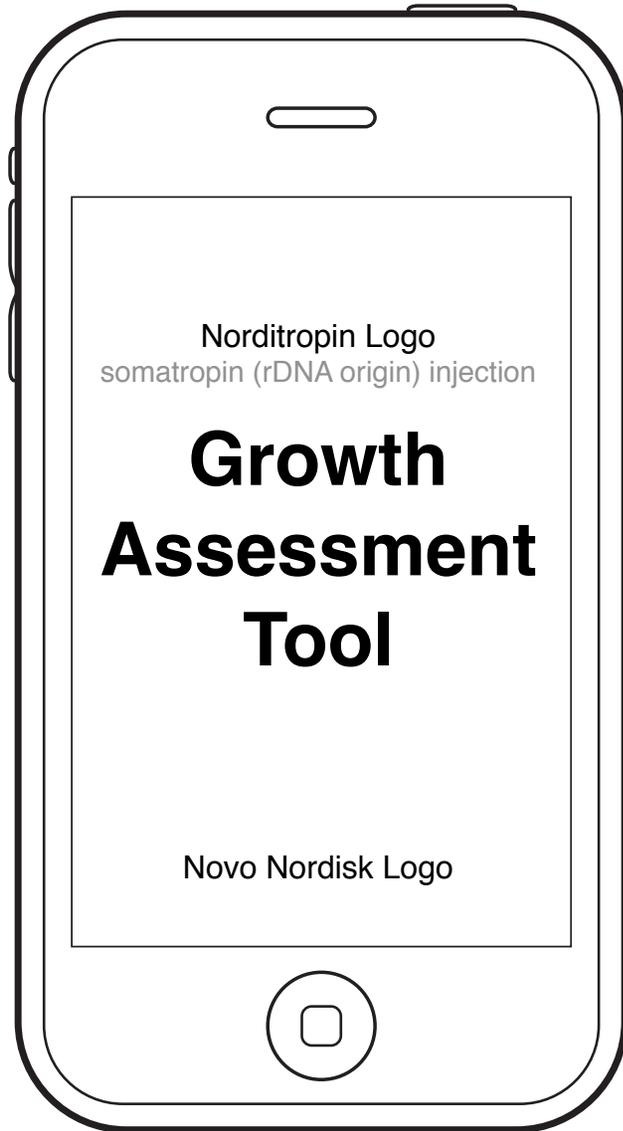


## Landing

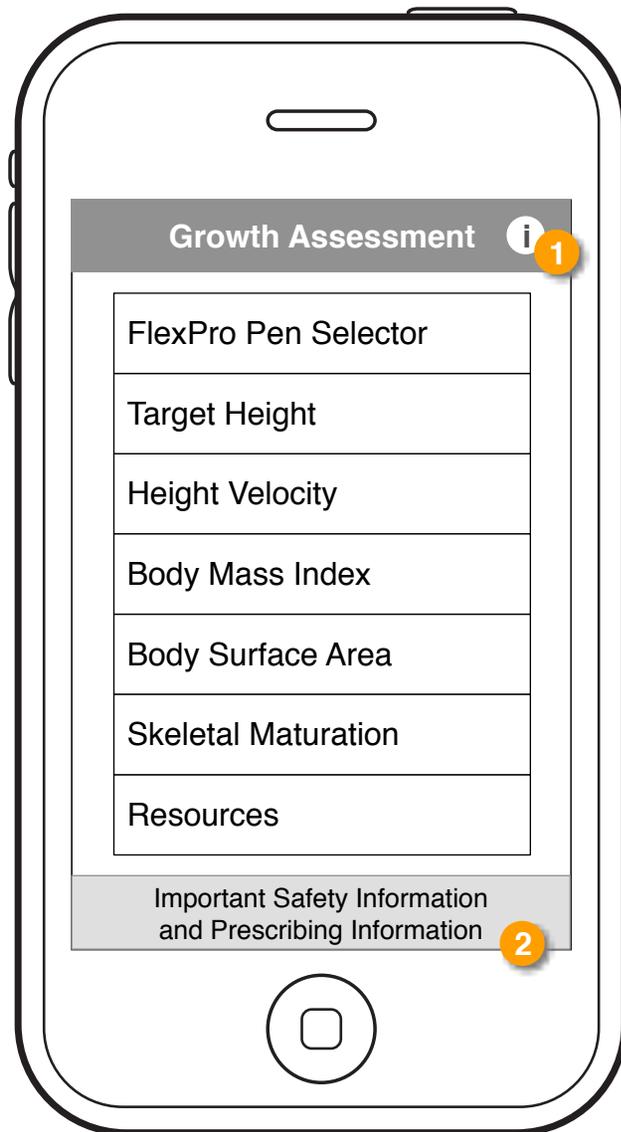


- 1** Information button: Tapping the "i" button opens the About screen (page 4).
- 2** Important Safety Information and Prescribing Information: Tapping the bar causes the ISI to slide up from the bottom (page 3). This element is global and always locked at the bottom of the screen, even when the page is scrolled up or down.

## Important Safety Information

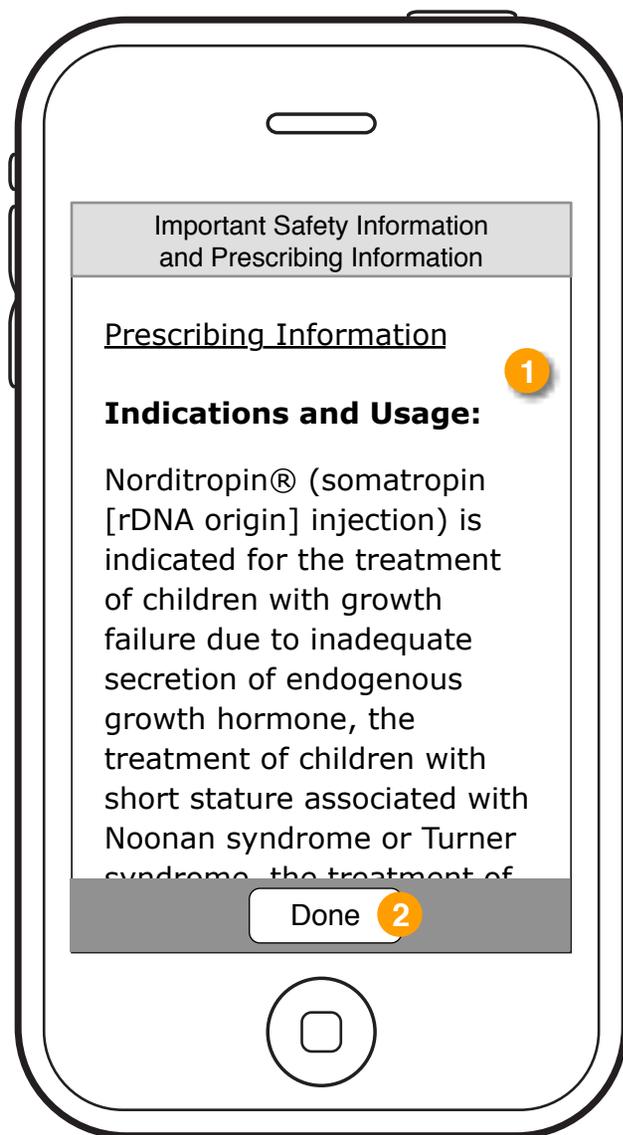


- 1 ISI: Users can scroll through the ISI. A link to the Prescribing Information can be found at the top.
- 2 Done button: Tapping Done closes the ISI and returns to the previous screen.

## About

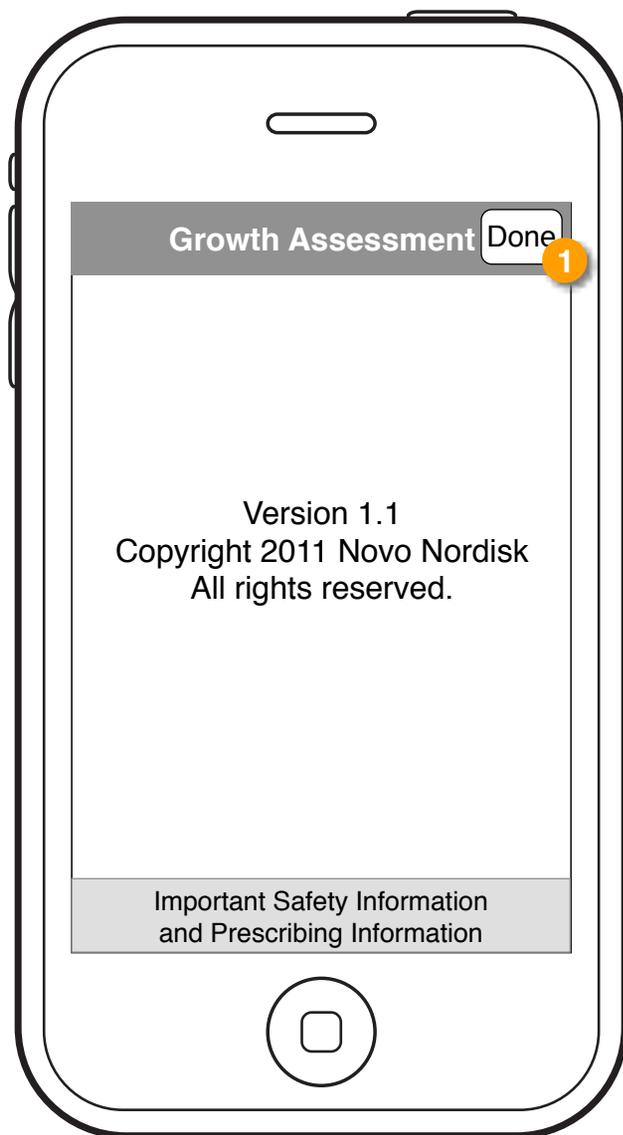

- 1 Done button: Tapping Done closes the About screen and returns to the landing screen (page 2).

## FlexPro Pen Selector

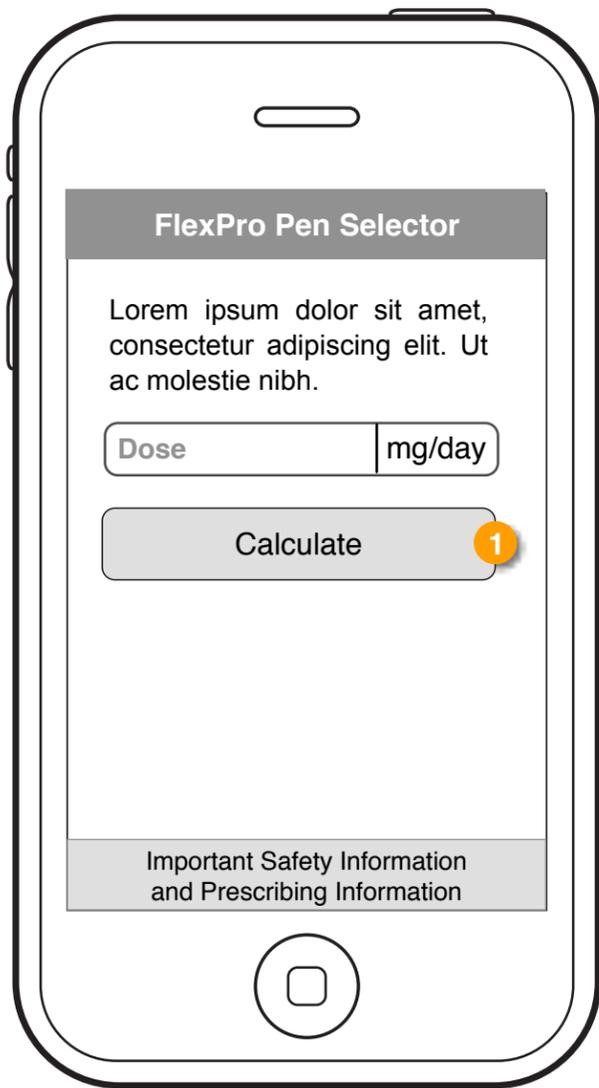

- 1 **Calculate:** Entering form data and tapping calculate will display the Results screen shown below.
- 2 **29-day supply:** This text will anchor link to the storage flexibility disclaimer at the bottom of the screen.
- 3 **Alert popup:** If the user enters a dose that cannot be injected by the pen, they will be given the option to round up or down to the nearest correct dosing increment.

FlexPro Pen Selector

Based on a dose of 1.5 mg/day:

image of pen

**The 5mg/1.5mL  
FlexPro pen will last  
3.3 days**

You will need 9.0 pens  
for a 30-day supply.\* 2

image of pen

**The 10mg/1.5mL  
FlexPro pen will last  
6.7 days**

You will need 4.5 pens  
for a 30-day supply.\* 2

image of pen

**The 15mg/1.5mL  
FlexPro pen will last  
10 days**

You will need 3 pens  
for a 30-day supply.\* 2

## Results

//

\* All Norditropin® products must be refrigerated prior to first use. Do not freeze. After initial use, FlexPro® 5 mg/1.5 mL and 10 mg/1.5 mL delivery pens can either be stored outside of the refrigerator (at up to 77°F) for use within 3 weeks, or in the refrigerator (between 36°F and 46°F) for use within 4 weeks. FlexPro® 15 mg/1.5 mL and Norditropin® NordiFlex® 15 mg/1.5 mL and 30 mg/3 mL delivery pens must always be refrigerated (between 36°F and 46°F)—both prior to and after the initial injection—for use within 4 weeks.

Important Safety Information and Prescribing Information

FlexPro Pen Selector

3

Note: You have entered a dose that needs to be rounded up or down to provide accurate results:

You entered .04 mg/day

Round up to .05 mg/day

Round down to .025 mg/day

Cancel

Important Safety Information and Prescribing Information

## Target Height

Target Height

Lorem ipsum dolor sit amet, consectetur adipiscing elit. Ut ac molestie nibh.

Mother's height	inches <sup>1</sup>
Father's height	inches

Male    Female

Calculate <sup>2</sup>

Important Safety Information and Prescribing Information

- 1 inches:** Tapping the inches button will toggle the units to metric (cm).
- 2 Calculate:** Entering form data and tapping calculate will display the Results screen shown below.
- 3 Units of measurement:** The results will display both standard and metric units. The top/larger result will depend on which unit was used when entering the information. For example, in this instance the user entered the height in inches, so the larger result is in inches with the smaller result below in cm.

## Results

Target Height

Mother's height: 66"  
Father's height: 73"  
Gender: Male

TARGET HEIGHT

**71"** <sup>3</sup>

67"    (+/- 4 inches)    75"

135 cm  
(125 - 145 cm)

Important Safety Information and Prescribing Information

## Height Velocity

Height Velocity

Lorem ipsum dolor sit amet, consectetur adipiscing elit.

Previous visit

Date

Height inches

Current visit

Date January 1, 2011 **1**

Height inches **2**

Calculate **3**

Important Safety Information and Prescribing Information

- 1** Date for Current Visit: This date field will default to today's date.
- 2** inches: Tapping the inches button will toggle the units to metric (cm).
- 3** Calculate: Entering form data and tapping calculate will display the Results screen shown below.
- 4** Alert popup: If the user enters a time period of less than 3 months, they will be alerted that results may not be accurate.
- 5** Continue: Tapping Continue will display the results screen with accuracy disclaimer.

## Results

Height Velocity

January 6, 2010: 62 inches  
August 24, 2010: 62.75 inches

HEIGHT VELOCITY

**1.25**  
inches/year

5 cm/year

Important Safety Information and Prescribing Information

Height Velocity

Warning: The time period you have entered is less than 3 months, and therefore may yield less accurate results. **4**

Cancel Continue

Important Safety Information and Prescribing Information

Height Velocity

January 6, 2010: 62 inches  
March 18, 2010: 62.75 inches

HEIGHT VELOCITY

**3.25** **5**  
inches/year

Warning: The time period you have entered is less than 3 months,

Important Safety Information and Prescribing Information

## Body Mass Index

Body Mass Index

Lorem ipsum dolor sit amet, consectetur adipiscing elit. Ut ac molestie nibh.

Weight	lbs
Height	inches

Calculate

Important Safety Information and Prescribing Information

**1** Units toggle: Tapping the lbs or inches button will toggle to metric. Tapping the years button will toggle to months.

**2** Calculate: Entering form data and tapping calculate will display the Results screen shown below.

## Results

Body Mass Index

Weight: 78 lbs  
Height: 56.5 inches  
Age: 6 Gender: Male

BODY MASS INDEX

**22.5**  
lb/in<sup>2</sup>

Important Safety Information and Prescribing Information

## Body Surface Area

Body Surface Area

Lorem ipsum dolor sit amet, consectetur adipiscing elit. Ut ac molestie nibh.

Weight	lbs
Height	inches

Calculate

Important Safety Information and Prescribing Information

**1** Units toggle: Tapping the lbs or inches button will toggle to metric. Tapping the years button will toggle to months.

**2** Calculate: Entering form data and tapping calculate will display the Results screen shown below.

## Results

Body Surface Area

Weight: 78 lbs  
Height: 56.5 inches  
Age: 6

BODY SURFACE AREA

**1.94**  
meters square

Average BSA for:  
adult men: 1.9 m2  
adult women: 1.6 m2  
9 year-old: 1.07 m2  
10 year-old: 1.14 m2  
12-13 year-old: 1.33 m2


Important Safety Information and Prescribing Information

## Skeletal Maturation




- 1 Male/Female: Male will be selected by default. User can toggle to Female.
- 2 Age selector: The user can scroll side to side to choose a specific age.

## Resources

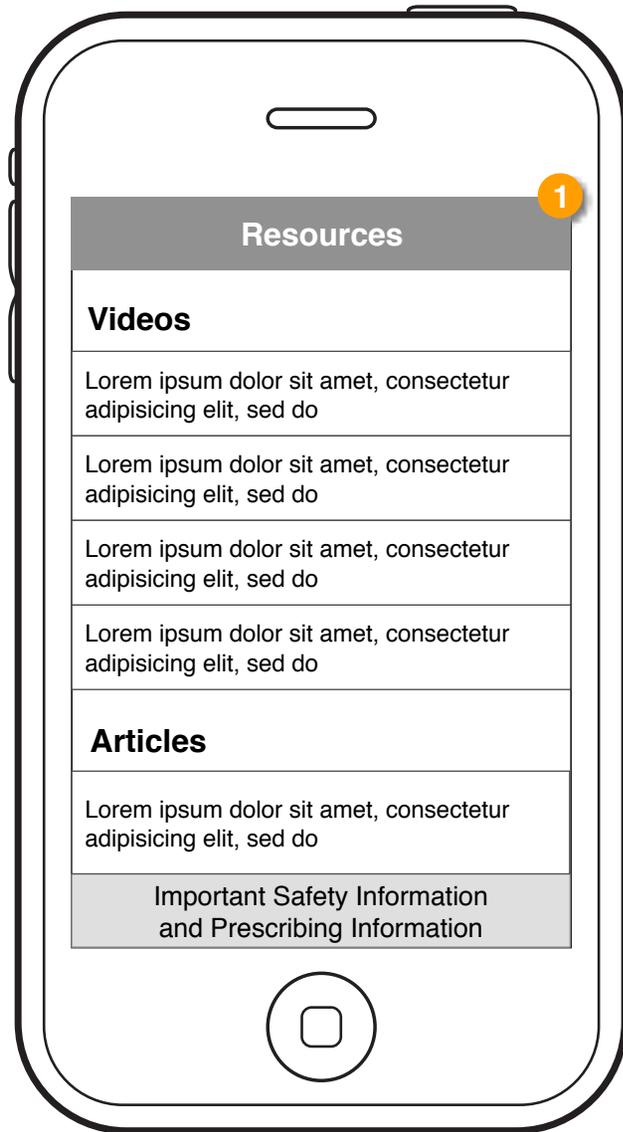

- 1 Resources: Displays Norditropin branded resources such as MOA, videos, vis-aids, etc.

## Input Methods

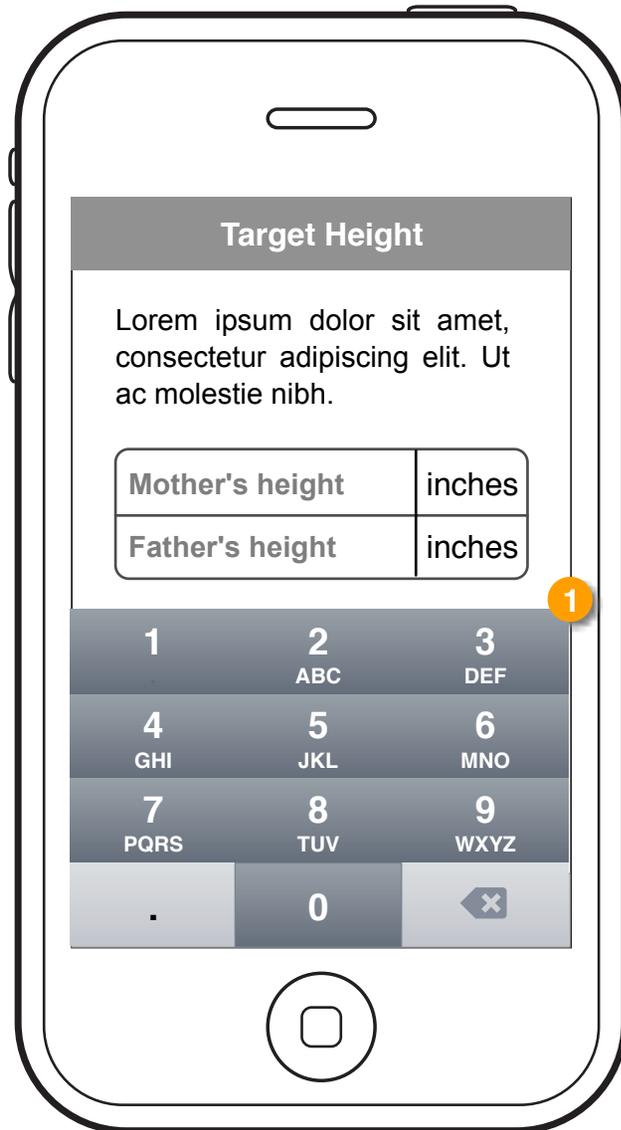

- 1 Data Input: The number pad will slide up from the bottom when a user taps a field to input numeric data, such as height, weight, etc. A date wheel will be used to enter any dates.

1

January	23	2008
February	24	2009
March	25	2010
April	26	2011
May	27	2012